




Refer to the user manual for comprehensive information:
<http://www.mttts-asia.com/wallaby-warmer/>


STEP 1 - Remove components from the box

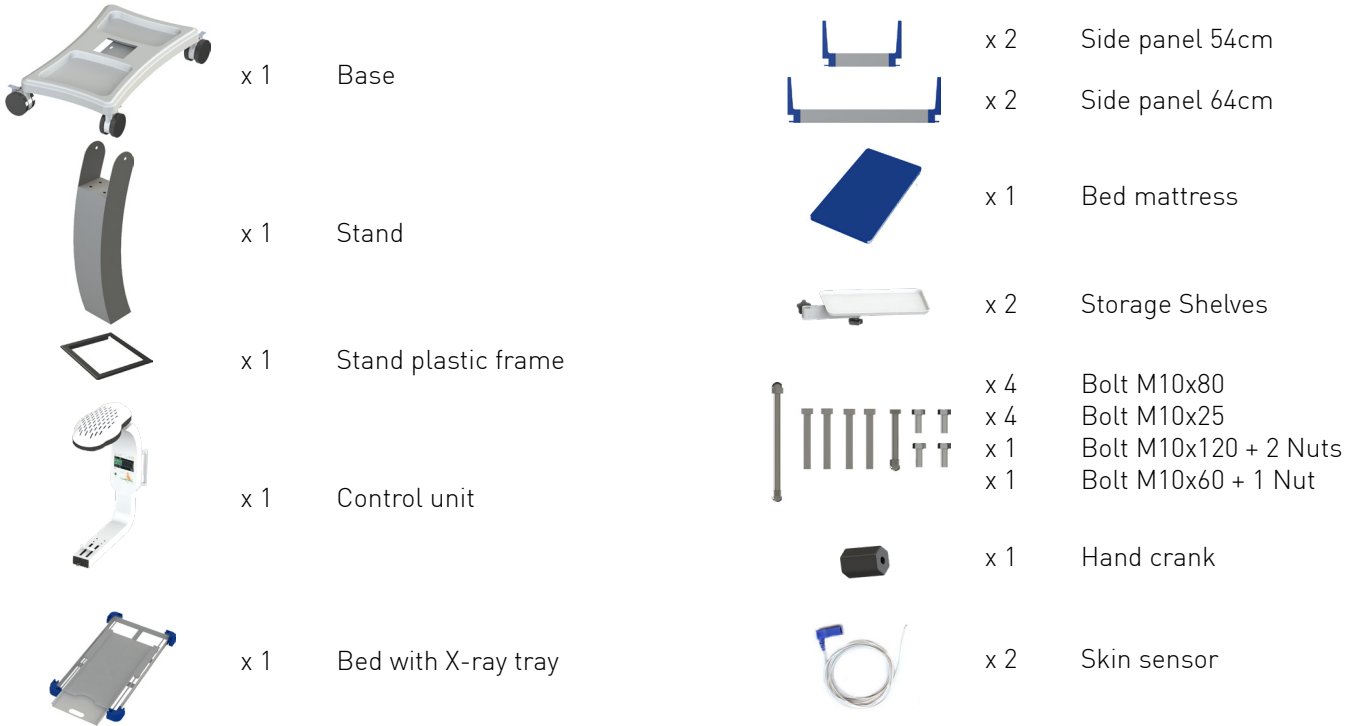


STEP 2 - Attach the stand to the base




Two people might be required to perform this step

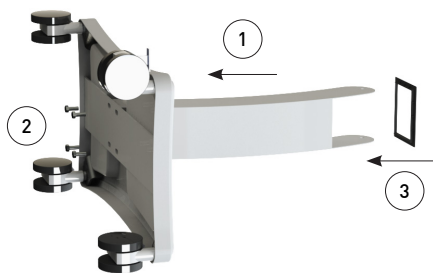

STEP 3 - Attach the control unit to and the hand crank to the stand

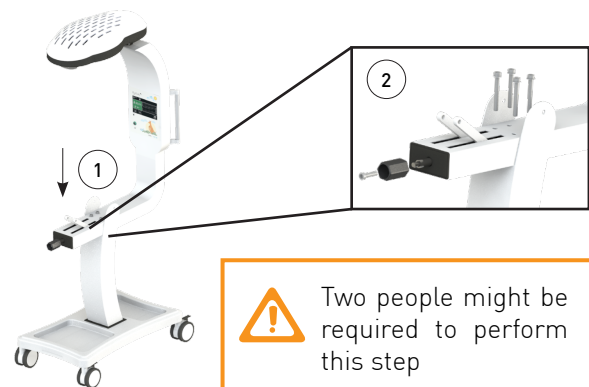


STEP 4 - Attach the side panel to the bed

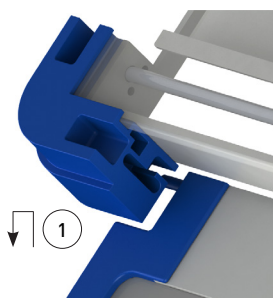


Attach only the side panel closest to the control unit in this step. Remaining panels will be attached in Step 6.

STEP 5 - Attach the bed

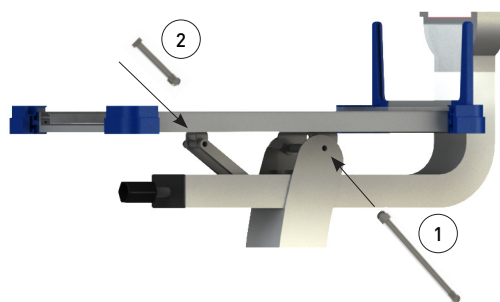


STEP 6 - Attach three remaining side panels to the bed

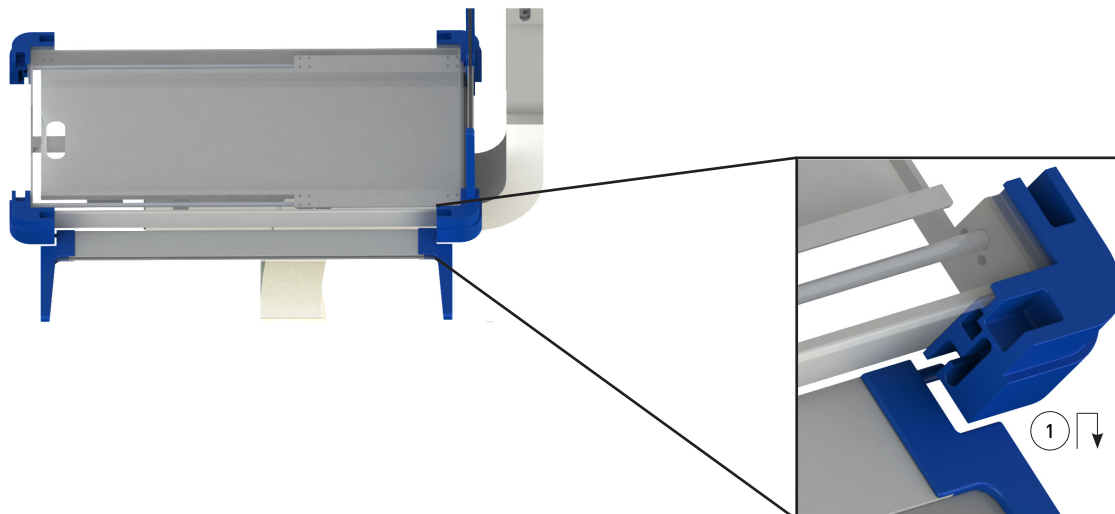

STEP 7 - Attach the storage shelves



STEP 8 - Plug in and power up

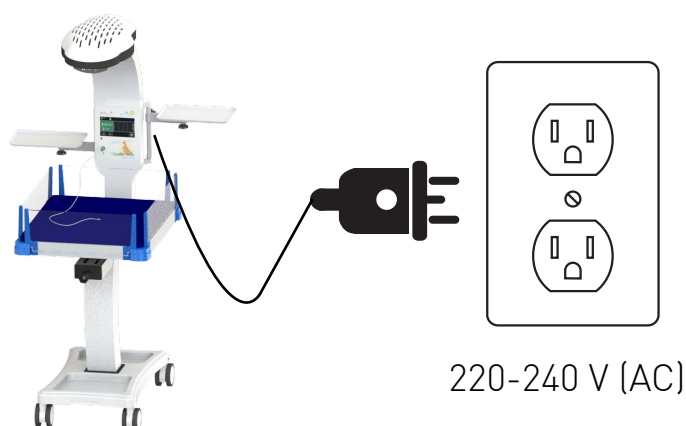

Make sure to remove the protector from the heater element before the first use

