



Refer to the user manual for comprehensive information:
<http://www.mttts-asia.com/wallaby-warmer/>

Optional components list



IV pole (ID 1331)



Storage drawer (ID 1332)



Cylinder holder (ID 1333)



Storage shelf (ID 1254)



Holding pole (ID 1792)


Cylinder holder (ID 1333)




2 x Bolt M10x25

Storage drawer (ID 1332)

- ① Disconnect top and bottom parts.
Press black latch on both sides and slide out.

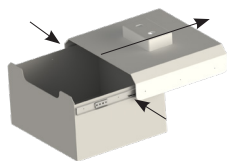

- ② Attach the top part to the control unit.
Use 2 x Bolt M10x25.

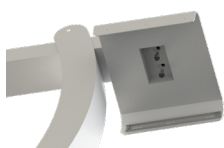

- ③ Connect the bottom part.
Slide it in until it clicks.



IV pole (ID 1331), storage shelf (ID 1254), holding pole (ID 1792)

