

Firefly 
phototherapy

LifeKit by **MTTS**



PHILIPS
LED

CE 2265

FIREFLY

A **double-sided phototherapy device** designed to treat **neonatal jaundice** with optional 2-hour battery back up.



Overall

- compact form fits in an infant cart or mother's bed
- handholds facilitate easy portability
- tight seams to prevent dirt build-up
- curved surfaces for easy cleaning
- sealed to keep out bugs, dust, and liquids
- no moving parts, no internal fans
- optional 2h battery back up

1. Top Light

- fixed at the most clinically effective height
- slim form for easy patient observation
- curved top prevents placement of objects
- soft edges to protect infant during placement

2. Bottom Light


- durable thick-walled plastic
- tilted bottom directs heat off back edge of device

3. Bassinet

- comfortable handhold around the rim
- removable for patient transport and diagnostics
- smooth surface for quick wipe-down
- clear scratch-resistant plastic for easy observation

4. Control Panel

- single power button controls both lights
- single light intensity setting
- treatment hours and total device hours display

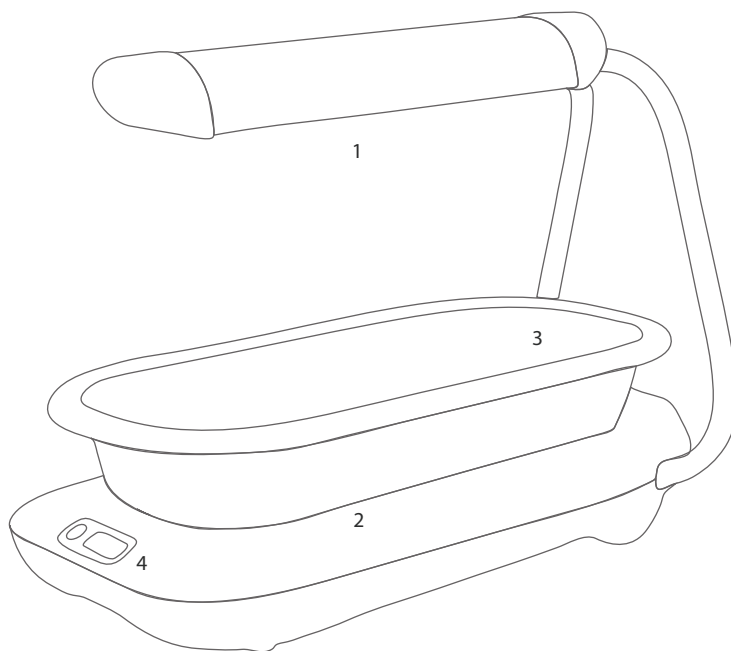

Dimensions and Weight

Unit dimensions (overall)	66 cm x 38 cm x 49.5 cm
Unit dimensions (top light)	51.6 cm x 12.5 cm x 6.3 cm
Unit dimensions (bottom light)	64.7 cm x 38 cm x 10 cm
Shipping (boxed) dimensions	57 cm x 43 cm x 70 cm
Net weight	13 kg
Shipping weight	17 kg

Performance Specifications

Type of LEDs	Philips LUXEON high-powered blue
Peak wavelength	465 to 485 nm
Lamp duration	60,000 hours
Average spectral irradiance	Top light - 40.0 $\mu\text{W}\cdot\text{cm}^{-2}\cdot\text{nm}^{-1}$ Bottom light - 66.0 $\mu\text{W}\cdot\text{cm}^{-2}\cdot\text{nm}^{-1}$
Effective surface area	Top light - 48 cm x 20 cm Bottom light - 48 cm x 20 cm
Irradiance uniformity ratio	Top light - 0.63 (IEC Compliant \rightarrow 0.4) Bottom light - 0.64 (IEC Compliant \rightarrow 0.4)

Electrical and Environmental Specifications

Power characteristics	60W, 100-240VAC, 47/63Hz
Off-the-shelf	ANSI/AAMI/EN 60601-1, UL/TUV
External Power Supply	Class I construction is standard (ground required) 100k hours MTBF Energy star efficiency level V compliant Over voltage / over current protection
Operating	Ambient temperature: +10°C to +35°C Humidity: 0% to 90% RH non-condensing Atmospheric pressure: 70-106kPa Sealing: IPX3
Transport and Storage	Ambient temperature: 0°C to +50°C Humidity: 0% to 90% RH non condensing Atmospheric pressure: 70-106kPa
Exclusion	Not for use inside an incubator/radiant warmer

Requirements

Agency Requirements	EN 60601-1:2006/AC2010	EN 60601-6:2010
	EN 60601-1-2:2007/AC2010	ISO 10993-Part 10
	EN 60601-2-50:2009	ISO 10993-Part 5

Ordering Information and Servicing

MAIN UNIT	ITEM ORDER NO	
Firefly Unit	FF1-1000	Main Unit
OPTIONAL ITEMS	ITEM ORDER NO	
Eye Patch	ID 564/5/6	Size: S, M, L (YPT-1010/1/2)
Battery Back up	ID 1494	Internal - 2h
Trolley	ID 1462	
Gel Mattress	ID 1525	
Lightmeter	LM-800	
User Manual		EN, VN, FR, SP, RU
WARRANTY		
Standard Warranty	1 year	Subject to MTTs policy



MEDICAL TECHNOLOGY TRANSFER AND SERVICES Co., LTD
ISO 13485: 2016 Certified Company

For more information or to place an order please contact us at:

info@mtts-asia.com

<http://mtts-asia.com>

+84 24 37 666 521

No 26, Lane 41 An Duong Vuong Street
Phu Thuong Ward, Tay Ho District
Hanoi, Vietnam

Distributor: